


MINIMAL STRESS: U/S ENERGY

Minimal Stress Technology

The patented Minimal Stress technology is the only system on the market that measures the actual movement of the tip during surgery and controls the power delivered to the handpiece in order to maintain the exact tip stroke set.

Minimal Stress guarantees the safest cataract extraction, optimizing energy at the lowest efficient level.

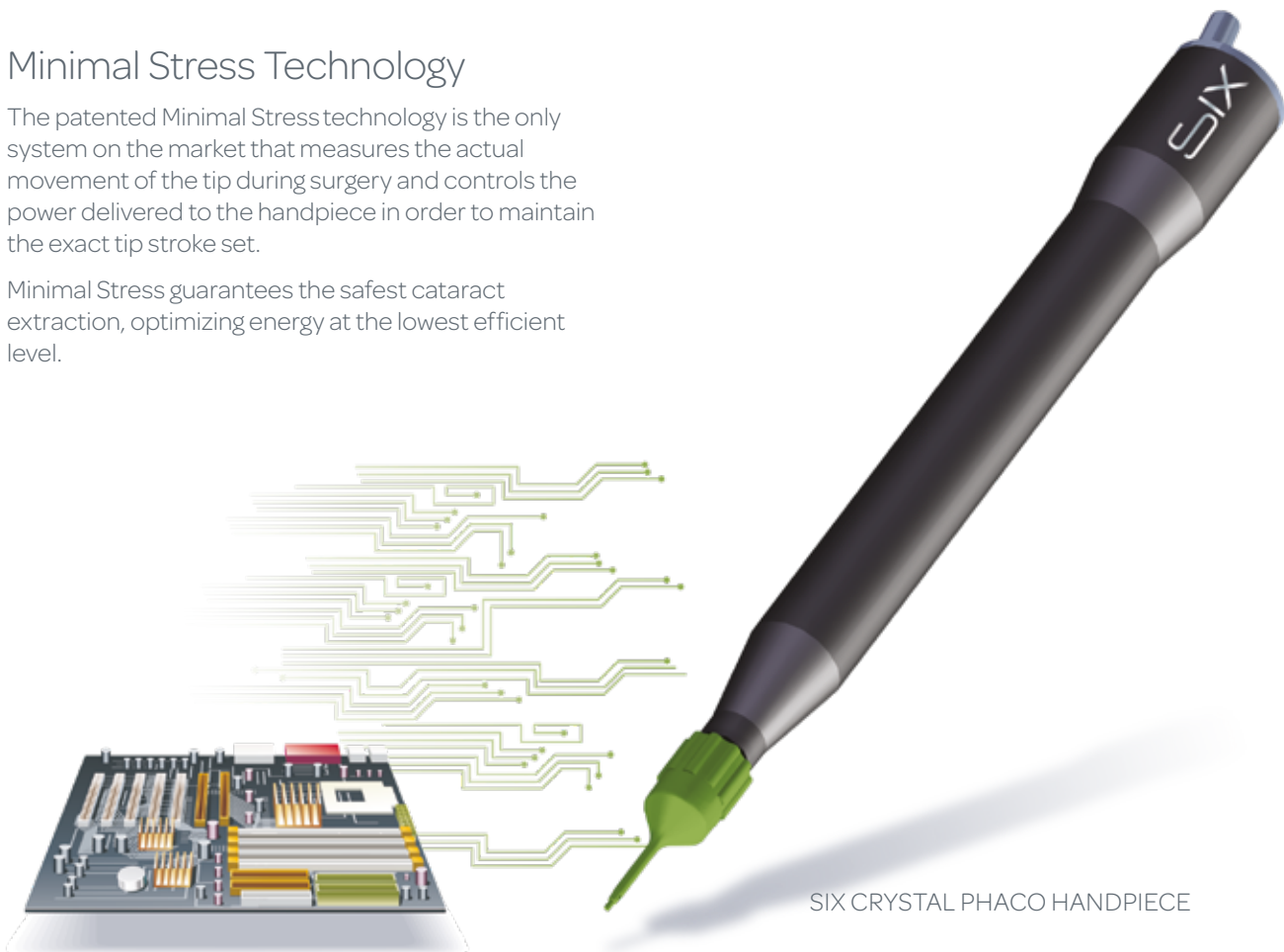

ELONGATION DYNAMIC WITH AND WITHOUT MINIMAL STRESS CONTROL

The chart shows the lower energy delivered into the eye by Minimal Stress control system compared to a standard U/S energy control system.

- Excess phaco energy without Minimal Stress control
- Elongation of the tip is not controlled
- Elongation of the tip is controlled by Minimal Stress

