

PREDICTIVE FLUIDICS




Double Pump System


The R-Evolution™ fluidics is managed by a double pump system, based on two latest generation pumps, one with flow control (peristaltic) and one with vacuum control (rotary vane for R-Evolution™ and Venturi for R-Evolution™ CR).

Easy choose between peristaltic and Venturi pumps through a single convenient I/A cassette.


The predictive fluidics ensures maximum flexibility and safety both in anterior and posterior segment surgery and with whichever technique used.

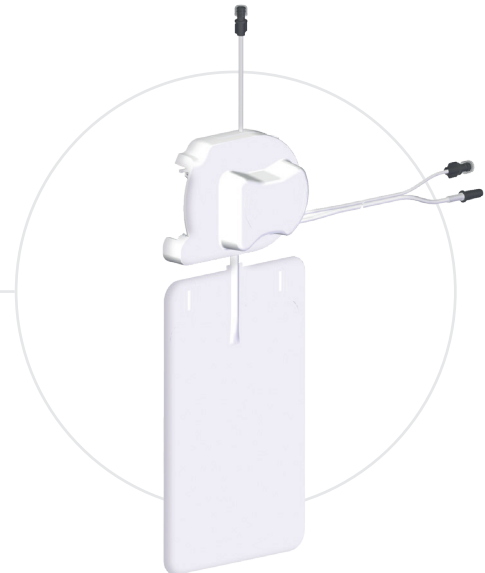

DOUBLE PUMP SINGLE PLATE WITH EASY ASSEMBLY I/A CASSETTE

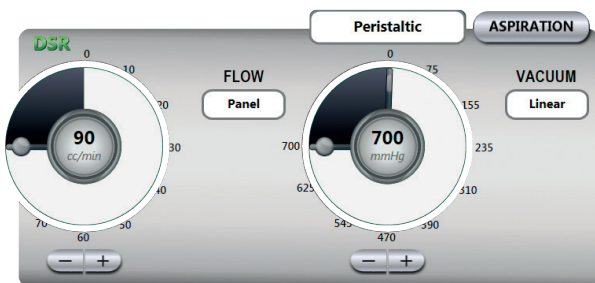

PERISTALTIC PUMP
Flow rate: up to 90 cc/min
Vacuum: up to 700 mmHg

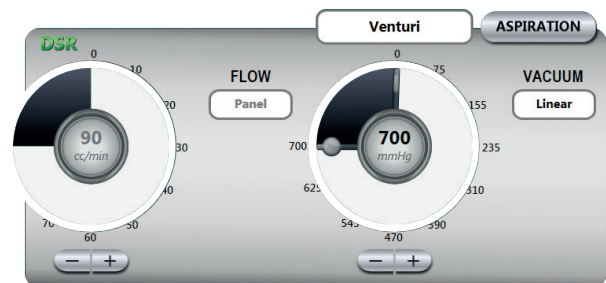

VENTURI/ROTARY VANE PUMP
Vacuum: up to 700 mmHg